

REVISION solutions

# GUARDIAN<sup>®</sup>

REVISION HINGE KNEE SYSTEM




**WRIGHT.**

Operative Flexibility • Component Modularity • Complete Instrumentation

# GUARDIAN<sup>®</sup>

REVISION HINGE KNEE SYSTEM

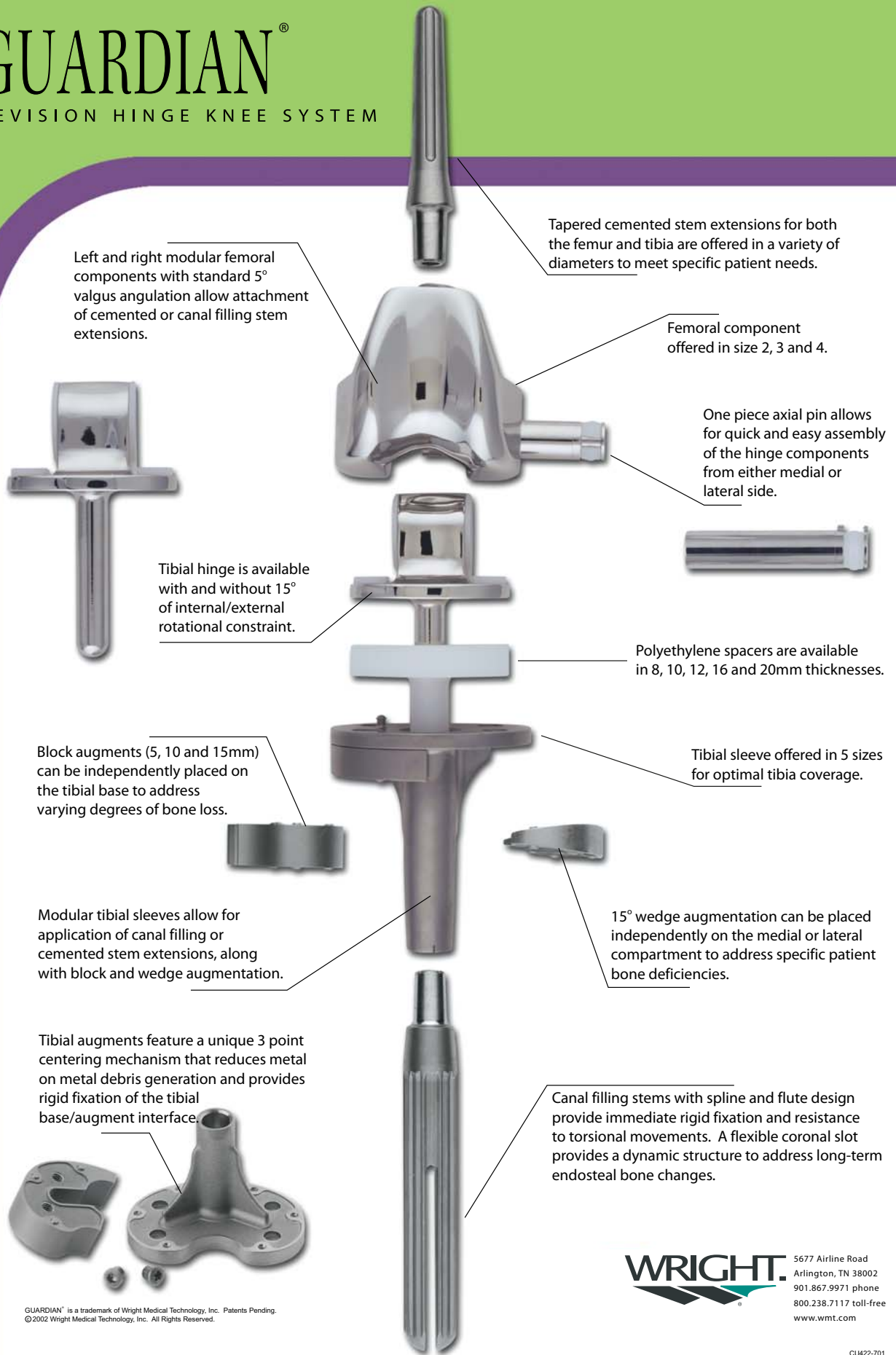

Left and right modular femoral components with standard 5° valgus angulation allow attachment of cemented or canal filling stem extensions.

Tapered cemented stem extensions for both the femur and tibia are offered in a variety of diameters to meet specific patient needs.

Femoral component offered in size 2, 3 and 4.

One piece axial pin allows for quick and easy assembly of the hinge components from either medial or lateral side.

Tibial hinge is available with and without 15° of internal/external rotational constraint.

Polyethylene spacers are available in 8, 10, 12, 16 and 20mm thicknesses.

Block augments (5, 10 and 15mm) can be independently placed on the tibial base to address varying degrees of bone loss.

Tibial sleeve offered in 5 sizes for optimal tibia coverage.

Modular tibial sleeves allow for application of canal filling or cemented stem extensions, along with block and wedge augmentation.

15° wedge augmentation can be placed independently on the medial or lateral compartment to address specific patient bone deficiencies.

Tibial augments feature a unique 3 point centering mechanism that reduces metal on metal debris generation and provides rigid fixation of the tibial base/augment interface.

Canal filling stems with spline and flute design provide immediate rigid fixation and resistance to torsional movements. A flexible coronal slot provides a dynamic structure to address long-term endosteal bone changes.

GUARDIAN<sup>®</sup> is a trademark of Wright Medical Technology, Inc. Patents Pending.  
© 2002 Wright Medical Technology, Inc. All Rights Reserved.

**WRIGHT.**  
5677 Airline Road  
Arlington, TN 38002  
901.867.9971 phone  
800.238.7117 toll-free  
www.wmt.com